

WorkiMed IVP Telescopic Light Cart



This new IVP medical trolley with telescopic pole is of light and robust aluminium construction and allows stable and intensive movement without strain.

Its large service basket 380 x 180 x 206 mm is fixed along the column at a height determined by the user.

The telescopic IVP pole made of anodised aluminium is adjustable in height from 700 mm to 1400 mm. The height is locked by a knurled ring.

Four hooks allow the suspension of medication bags.

The pole is fixed at the head of the column for a perfect adjustment of the height of the medication bags. Equipped with cable hooks to optimise the orderly storage of cables.

The movement is done with the carrying handle included in this configuration of WorkiMed IVP telescopic cart IVP Light Cart.

Equipped with comfortable CE medical anti-static wheels, the new medical cart can be manoeuvred effortlessly by the nursing staff at the point of care, the diameter of the wheels is 100 mm and they are all equipped with foot brakes. These wheels allow easy handling over bumps and imperfections in the floor.

* We can also make custom made trolleys according to your needs and in your hospital colours as an option









WorkiMed[®]
WORK PLACE TECHNOLOGY

Presentation




Part Number: WM 555.115

WorkiMed IVP Telescopic Light Cart

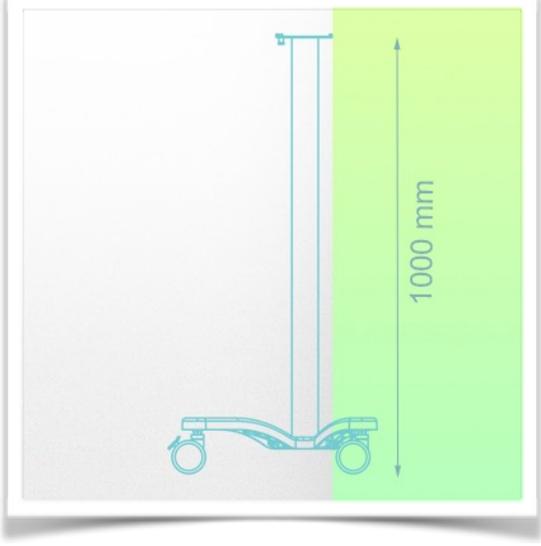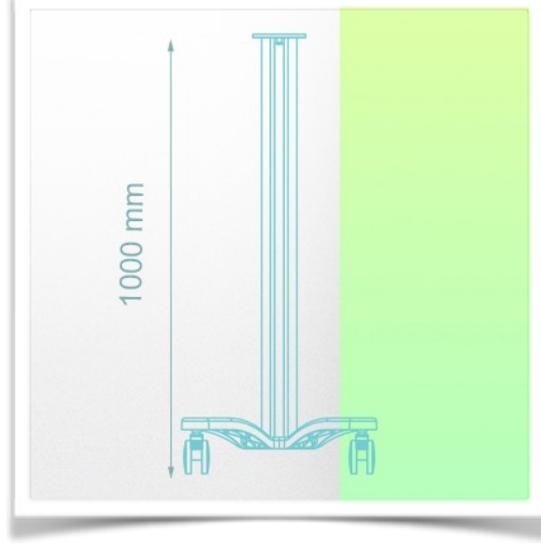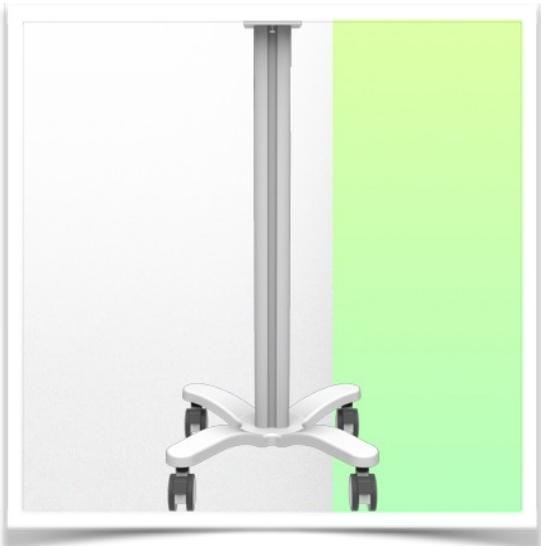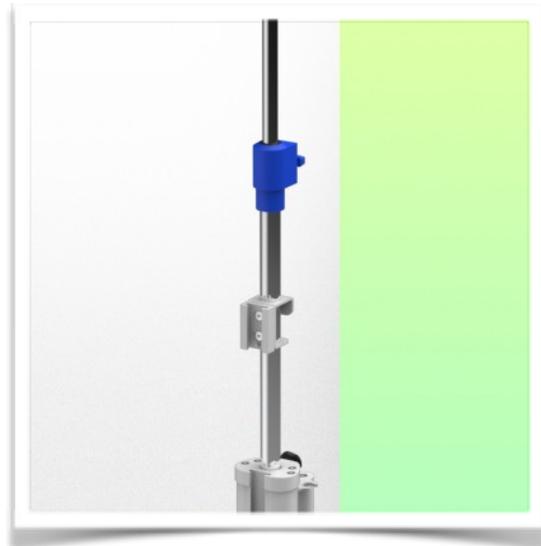









WorkiMed IVP Telescopic Light Cart

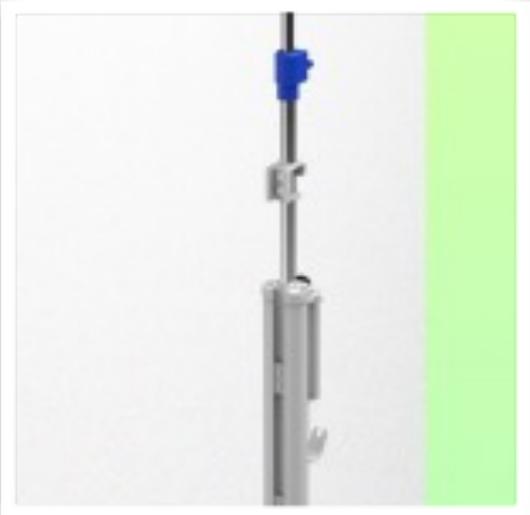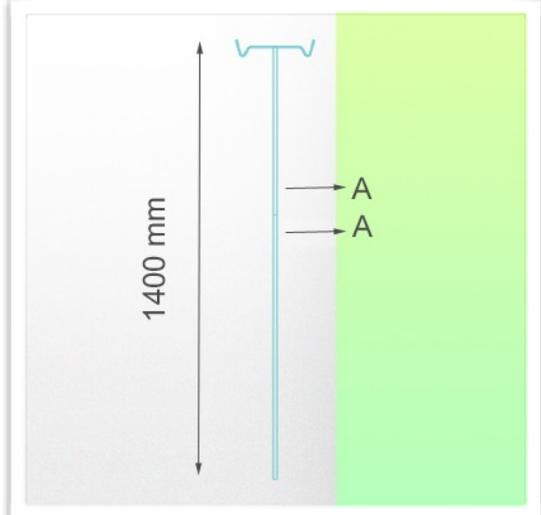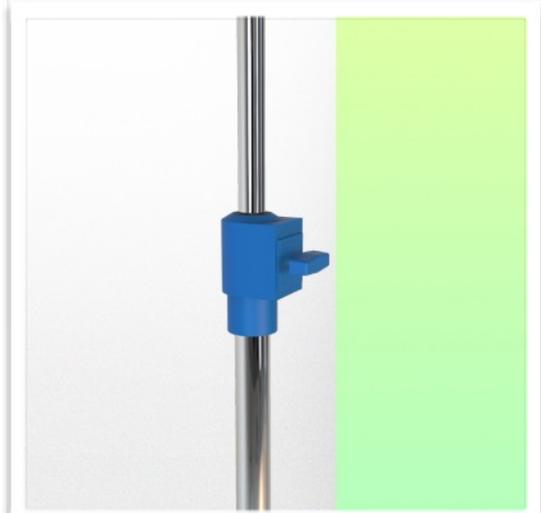




WorkiMed IVP Telescopic Light Cart

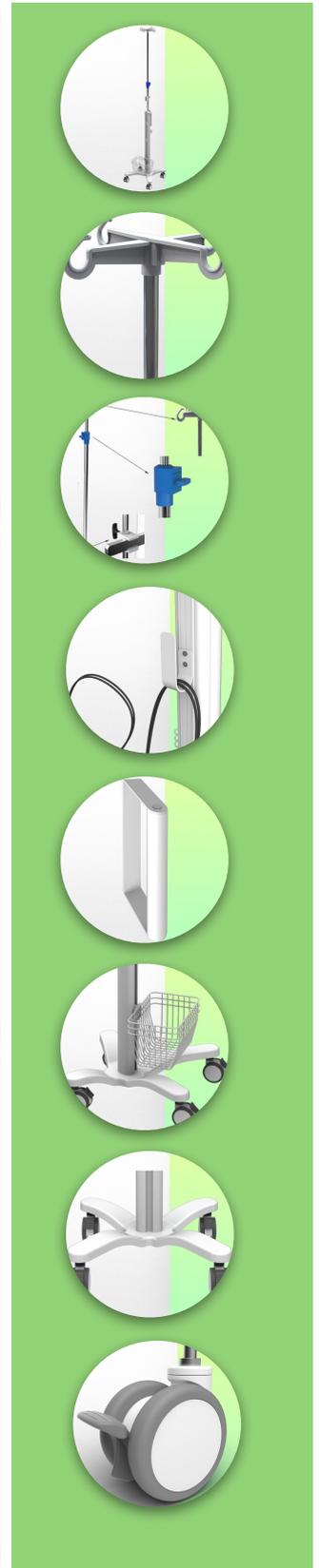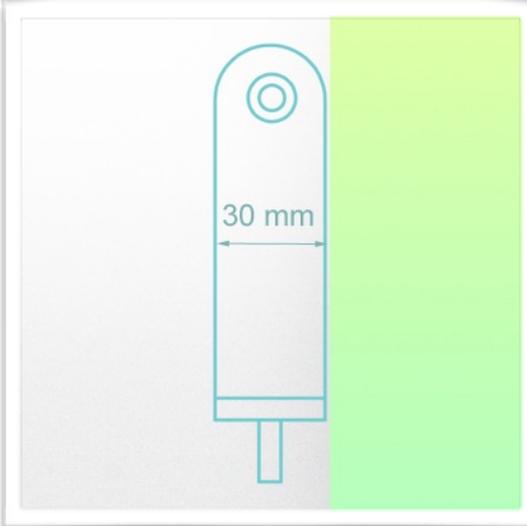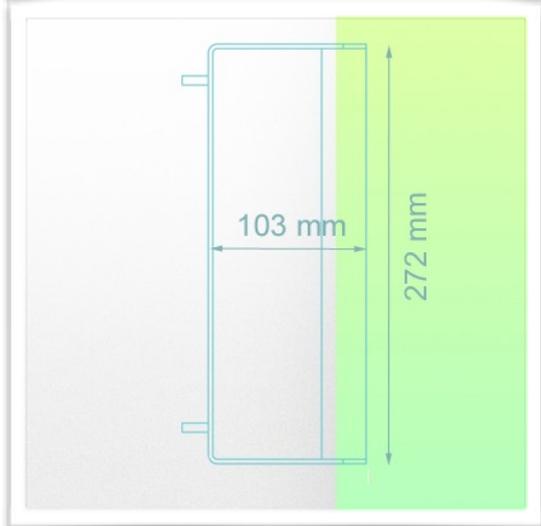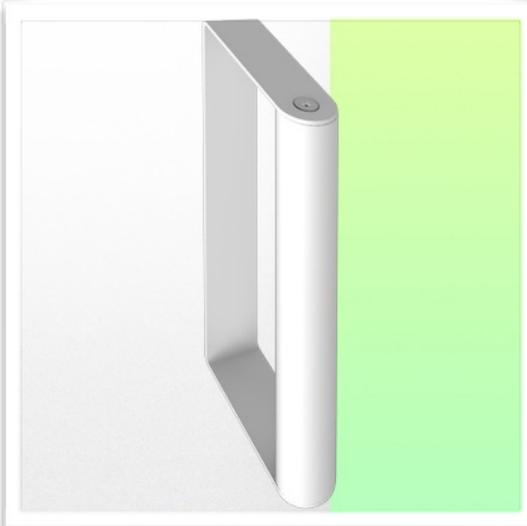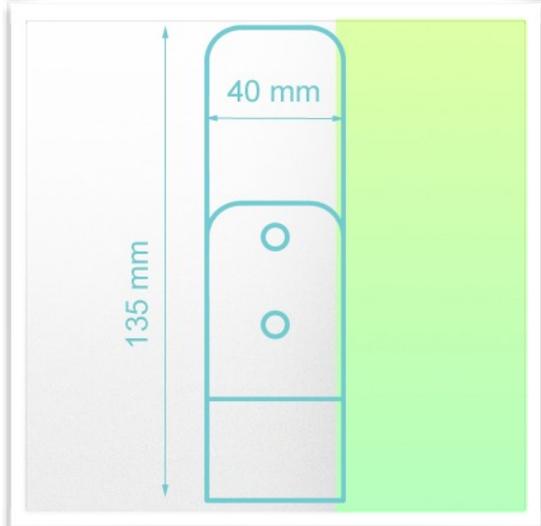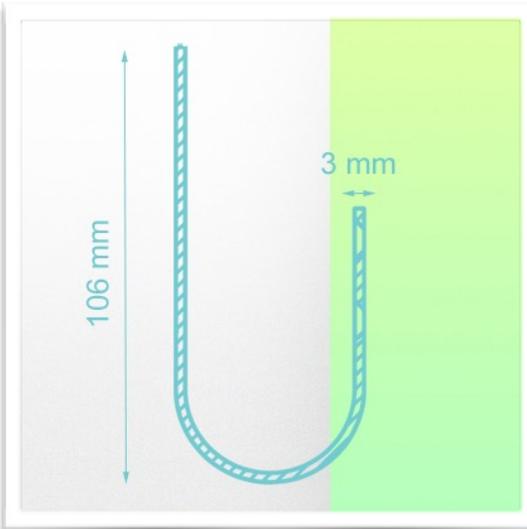



WorkiMed IVP Telescopic Light Cart

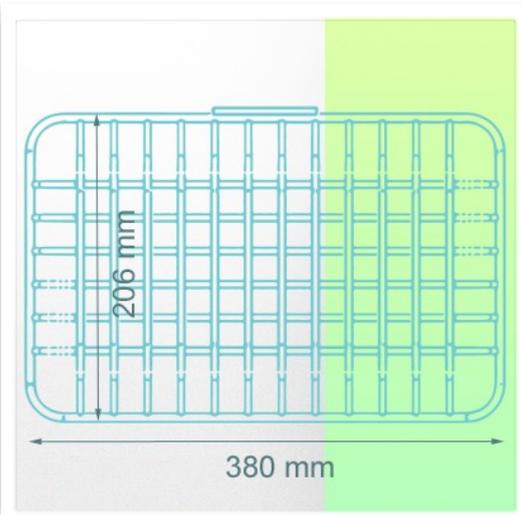




WorkiMed[®]
WORK PLACE TECHNOLOGY

Presentation

Part Number: WM 555.115

WorkiMed IVP Telescopic Light Cart

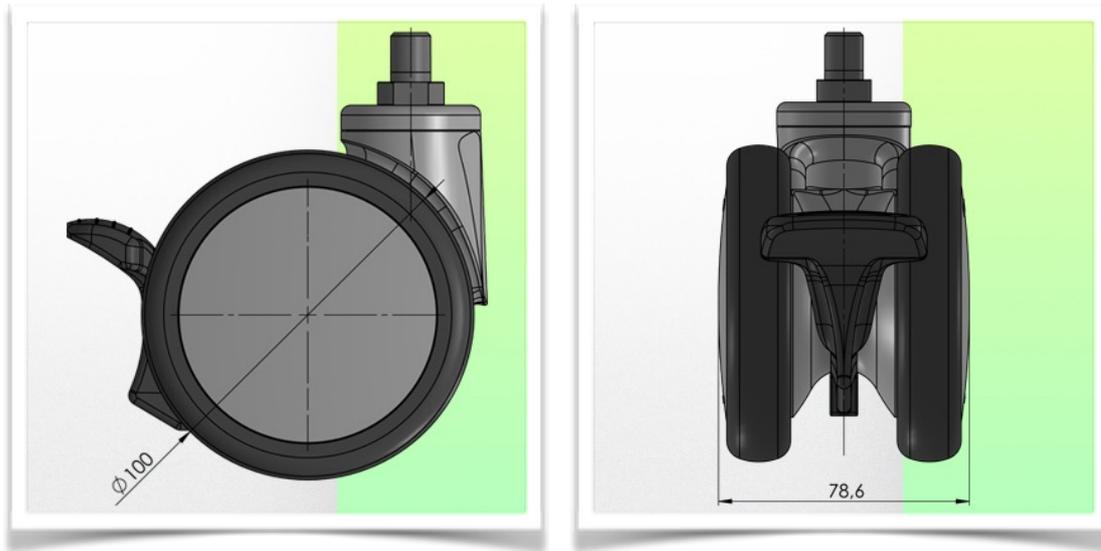

WorkiMed IVP Telescopic Light Cart

Technical specifications:

- * * Height of the column with base: 1000 mm
 - * Height adjustment capacity of the telescopic pole: 700 mm
 - * Wheel construction: Polyamide with polyurethane tread.
 - * Wheels/Lock: 4 locking wheels with brake, 100 mm diameter
 - * Trolley weight: 12 Kg
 - * Loading capacity: 20 Kg
 - * Main assembly: Aluminium
 - * Base diameter 630 mm
 - * Cable hook 106 mm
 - * Telescopic pole 700/700 mm
 - * 4 hooks for medicine bags
 - * Large basket 380 x 180 x 206 mm
 - * Integrated suspension hose management
 - * Carrying handle
 - * Colour: Ral 9016 White
 - * Certificates are in accordance with : IEC60601-1: 2012 (Edition 3)
 - * CE/Rohs compliant
 - * Fully recyclable
 - * 10 year warranty
 - * Wide range of accessories
 - * For faster and easier cleaning
 - * Simple and lightweight design
 - * Supplied fully assembled with minimal packaging
 - * Made in Europe (MCC John Gough design and development)
- * * We can also make custom-made trolleys according to your needs and in your hospital's colours as an option



WorkiMed IVP Telescopic Light Cart

Cleaning procedure

1. The trolley can be cleaned with most mild, non-abrasive solutions commonly used in hospitals (e.g. dilute bleach, ammonia or alcohol solutions).
2. The surface finish will be permanently damaged by strong chemicals and solvents such as acetone and trichloroethylene.
3. Do not use steel wool or other abrasive materials to clean the mounting assembly.
4. Damage caused by the use of unapproved substances or processes will not be covered by the warranty. We recommend that you test any cleaning solution on a small area that is not visible to check for compatibility.
5. Immediately wipe off any cleaning agent on the arm with a cloth dampened with water.
6. Dry the trolley thoroughly after cleaning.

CAUTION: WorkiWay SRL makes no claims regarding the effectiveness of the chemicals or processes listed as a means of infection control.

Consult your hospital infection control officer or epidemiologist. To clean or sterilize mounted instruments or accessories, refer to the specific instructions provided with those products.

Translated with www.DeepL.com/Translator (free version)

